

HD INTEGRATED ENDOSCOPY SYSTEM

STROBOLUX III LED+HD

- HD endoscopy camera
- Stroboscope
- LED light source



The Strobolux III LED + HD integrated endoscopy system is revolutionary in the ENT practice.

It allows the user to carry out stroboscopy and continuous light examination synchronized with the Optomic HD camera.

INTEGRATED SYSTEM STROBOLUX III LED+HD


HD ENDOSCOPY CAMERA

- HD output, HDMI type.
- 1080p HD digital resolution.
- Global shutter.
- Automatic white balance.
- Automatic synchronized ELC.
- Synchronization of camera and stroboscope to optimize the image quality.
- Camera head: 2 buttons for freeze and white balance.



STROBOSCOPE

- LED stroboscopic system.
- 2 independent and simultaneous light outputs.
- Phase and slow movement function.
- Analysis frequency from 60 to 1,000Hz.
- Phase range variation from 0° to 400°.
- Voice analysis: spectrogram and waveform.
- Air and contact microphone.
- 5" TFT display.

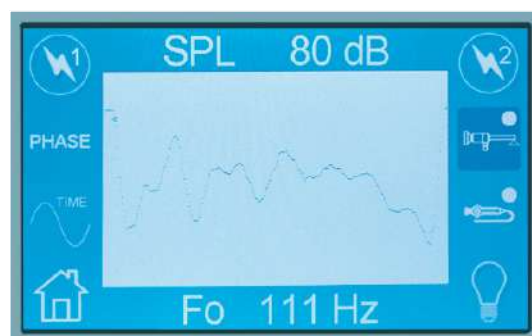

LED COLD LIGHT SOURCE

- LED type illumination.
- 2 independent outputs.
- Alternating or simultaneous function of the outputs.
- Life of LED lamp: 60,000h.
- Intensity of use memory.
- Fiber connection Storz type. Olympus type optional.



OPTOMIC ESPAÑA S.A.

Madroño, 4. E-28770 Colmenar Viejo. Madrid, Spain
Tel: +34 91 804 44 44
export@optomic.com
www.optomic.com



Specifications are subject to change without notice